



# Product Catalogue

2025



AVIAN  
BOVINE  
CANINE  
EQUINE  
FELINE  
PORCINE



Holsim

ANIMAL CARE & EDUCATION

# CONTENT PAGE

---

01 —

## AVIAN

|                               |   |
|-------------------------------|---|
| Kea Bandaging Model           | 4 |
| Cockatoo Bandaging Model      | 5 |
| Scarlet Macaw Bandaging Model | 6 |

02 —

## BOVINE

|                                |    |
|--------------------------------|----|
| Interlocking Bovine            | 8  |
| Birthing Cow Fibreglass        | 9  |
| Birthing Calf's                | 10 |
| Birthing Shroud & Rump         | 11 |
| Uterus Bag                     | 12 |
| Inflatable Bag                 | 13 |
| Interlocking Bovine Add on Kit | 14 |
| Interlocking Bovine Add On Kit | 15 |
| Calf Disbudding Kit            | 16 |
| Calf Disbudding Horn Buds      | 17 |
| Calf Disbudding Tool           | 17 |
| Metri-Treating Model           | 18 |
| Milking Udder                  | 19 |
| Milking Teats Kit              | 20 |
| Mastitis Teats Kit             | 21 |
| Udder Anatomy Model            | 22 |
| Cow Tail with Blood Taking Kit | 23 |

03 —

## CANINE

|                                      |    |
|--------------------------------------|----|
| Urinary Catheter Model (Female)      | 26 |
| Brachycephalic Intubation Model      | 27 |
| Canine Intubation Model              | 28 |
| Oesophageal Feeding Tube Model       | 29 |
| Canine Oesophageal Jugular Block     | 30 |
| Otoscope Ear Model                   | 31 |
| Jugular Blood Sampling Kit           | 32 |
| Cardio Audio Box & Lung Auscultation | 33 |
| IV Dog Leg Standard Kit              | 34 |
| IV Insert - Standard                 | 35 |
| Hind Limb IV Bandaging Leg Kit       | 36 |
| IV Insert - Bandaging Hind Limb      | 37 |
| Forelimb Canine Bandaging & IV Kit   | 38 |
| IV Insert - Limb                     | 39 |
| Articulated Hind Leg Bones           | 40 |
| Dog Bones                            | 41 |
| Mesocephalic Dental Dog              | 42 |
| Mesocephalic Dental Dog Teeth Set    | 43 |

04 —

## EQUINE

|                               |    |
|-------------------------------|----|
| Distal Limb Nerve Block Model | 45 |
| Nasogastric Intubation Model  | 46 |

05 —

## FELINE

|                                    |    |
|------------------------------------|----|
| Male Urinary Catheter Model        | 48 |
| Male Urinary Catheter Suture Block | 49 |

06 —

## PORCINE

|                                           |    |
|-------------------------------------------|----|
| Interlocking Porcine                      | 51 |
| Intramuscular Block with Blood Taking Kit | 52 |
| Pig Ear                                   | 53 |
| Ear Blood Taking Kit                      | 54 |
| Nose Ringing/Snaring Model                | 55 |
| Jugular Blood Taking Kit                  | 56 |
| AI Interlocking Porcine                   | 57 |

07 —

## ACCESSORIES

|                        |    |
|------------------------|----|
| Suture Pad             | 59 |
| Suture Pad Base        | 60 |
| Suture Kit             | 61 |
| Amniotic Fluid         | 62 |
| Sim Poo                | 63 |
| Blood Pump             | 64 |
| Artificial Blood       | 65 |
| Blower & Blower Tube   | 66 |
| Bungee & Carabiner Kit | 67 |
| Holsim Wooden Mallet   | 68 |

---



02 —

# Avian

*A survey showed that 95% of students felt that using training manikins had improved their psychomotor skills. This emphasises that when alternatives are used instead of live animals for training vet students, the same academic outcome can be reached or even exceeded.*

# Avian

## Kea Bandaging Model

**SKU**  
**A100-000**

**The Kea Bird Bandaging Model is an anatomically correct model for practising bandaging techniques in order to reduce pain and stabilise fractures.**

The Kea has humerus, radius, ulna, metacarpus and phalanx bones with primary, secondary and tertiary feathers attached.

Stretchy silicone connects the wing bones to the body to replicate the tendon's instinctive drawback. Positionable wings allow them to be rotated and fold against the body naturally.

The model is weighted accordingly to represent a juvenile male of the species selected.

*Other colours/bird breeds of similar size, as below, may be requested, or email us with the species you would like to see!*

---

**Replacement Feet Available**

---



# Avian

## Cockatoo Bandaging Model

**SKU**  
**A200-000**

**The Cockatoo Bandaging Model is an anatomically correct model for practicing bandaging techniques in order to reduce pain and stabilise fractures.**

The cockatoo has humerus, radius, ulna, metacarpus and phalanx bones with primary, secondary and tertiary feathers attached.

Stretchy silicone connects the wing bones to the body to replicate the tendon's instinctive drawback. Positionable wings allow them to be rotated and fold against the body naturally.

The model is weighted accordingly to represent a large male sulphur- crested cockatoo. The head can also swivel 360 degrees.

*Other colours/bird breeds of similar size, as below, may be requested, or email us with the species you would like to see!*

---

**Replacement Feet Available**

---



# Avian

## Scarlet Macaw Bandaging Model

**SKU**  
**A300-000**

**The Scarlet Macaw Bandaging Model is an anatomically correct model for practicing bandaging techniques in order to reduce pain and stabilise fractures.**

The Macaw has humerus, radius, ulna, metacarpus and phalanx bones with primary, secondary and tertiary feathers attached.

Stretchy silicone connects the wing bones to the body to replicate the tendon's instinctive drawback. Positionable wings allow them to be rotated and fold against the body naturally.

The model is weighted accordingly to represent a large male scarlet macaw. The head can also swivel 360 degrees.

*Other colours/bird breeds of similar size, as below, may be requested, or email us with the species you would like to see!*

---

**Replacement Feet Available**

---





02 —

# Bovine

*Animal training models, simulators and manikins allow for repetitive practice, even of procedures that are clinically uncommon or associated with high risk, in environments that are comparatively 'low stakes' and free from stress. They also offer standardisation of the learning experience and assessment processes for students.*

# Bovine

## Interlocking Bovine

**SKU**  
**B100-000B**

**Holsim's Interlocking Bovine is a flat-pack, transportable training model for classroom teaching of calving techniques. This interlocking birthing cow is made from plywood pieces that slot and bolt together for easy setup and pack down.**

Veterinary students are provided with the opportunity to experience the calving process in a controlled environment with clear assessment and feedback from the tutor. Once the student has this initial training, they are well prepared to get the most out of live animal training and the experiences that follow.

The top cage of the bovine is easily lifted off to set up and pack down the uterus bag, calf, etc.

Both the Holsim Milking Udder and Cupping Udder are compatible with this model and can be bolted onto the underside of the bovine for practising milking, mastitis treatment applications and cupping.

The interlocking bovine has carry handles on both the hinged base and the front of the bovine for moving it easily.





# Bovine

## Birthing Cow — Fibreglass

**SKU**  
**B1100-000**

**A transportable training model for classroom teaching of calving techniques. This life-size bovine enables students to practise manipulating and birthing the calf safely in a controlled environment.**

An inflatable bag acts as soft tissue and muscle to support the uterus bag in which the calf is zipped inside of with our Holsim-Amniotic Fluid. This enables lubrication during the birthing procedure. Modelled from actual bovine anatomy, each cow has its own uterus, rectum and birth canal.

The Holsim Birthing Bovine is heated internally, and has a blow-back switch for the rectum to replicate gas build up, enabling students to learn how to relieve this. For ease of cleaning and maintenance, the uterus bag unzips from the birth canal, whilst the rump and tail can be removed from the cow. Clear windows in the sides and top of the bovine allow teachers to see how students are going, and give live feedback during training.



## Bovine Birthing Calf

### SKU

B100-200J - Jersey

B100-200A - Angus

B100-200F - Friesian

**The Holsim Birthing Calf is a transportable training manikin to be used with the Holsim Interlocking Bovine and fibreglass, for classroom teaching of calving techniques.**

Modelled from actual bovine calf anatomy, the life-like calf can be placed in the Interlocking Bovine's uterus in a range of positions. This allows students to experience both easy and complicated procedures where the calf needs to be manipulated in the uterus before birthing.

The calf has internal bones in the hips and legs providing accurate movement and skeletal limits of the legs. This ensures a realistic challenge to the student when manipulating the position of the calf and the strength to apply force during the birthing procedure.



# Bovine

## Birth Shroud & Rump

**SKU**  
**B100-300**

**For use with the Holsim Interlocking Bovine and Holsim Fibreglass Birthing Cow.**

This silicone rectum, vaginal entrance and rump are attached to a hard shroud that mounts onto the Holsim Interlocking Bovine backside using threaded knurl screws to secure it in place.

The large vaginal opening allows for the Holsim Birthing Calf to be pulled through easily. The Uterus Bag zips onto the vaginal entrance for a smooth transition between the two parts. The Holsim Rectum Bag can attach to the rectum ring.



# Bovine

## Uterus Bag

**SKU**  
**B100-400**

**Designed to be used with the Holsim Interlocking Bovine for birthing simulation and Holsim Birthing Cow.**

This is a durable, waterproof bag designed to be used in birthing simulation scenarios with the Interlocking Bovine.

The bag has a zip which allows you to insert the Holsim Birthing Calf easily. Holsim Amniotic Fluid can be applied inside the bag over the calf.

The bag has suspension holes for attaching bungee cords/carabiners to suspend the bag inside the cow.



## Bovine Inflatable Bag

**SKU**  
**B100-500**

**Designed to be used with the Holsim Interlocking Bovine for birthing simulation.**

This is a durable, inflatable bag designed to support the calf in birthing simulation scenarios.

The bag can be inflated with the Holsim Blower. This bag supports the Uterus Bag and Calf from underneath, providing a soft cushioned surrounding for the calf to be enclosed in.



1a —

## INTERLOCKING BOVINE

### Bovine Birthing Kit

#### SKU B100-K100

- Birthing Calf with Bones
- Uterus Bag
- Inflatable Bag
- Blower & Blower Tube
- Birthing Shroud Mount
- Birthing Shroud Rump & Vaginal Entrance
- Rectum Bag & Clamp
- Tail
- Amniotic Fluid
- Bungee & Carabiner Kit
- Pelvis Cushion

When practising the birthing process, the Birthing Calf should be covered in Amniotic Fluid while being placed in the uterus bag. This represents the amniotic sac as realistically as possible and allows the calf to slide out of the birth canal naturally and smoothly.

### Artificial Insemination Kit

#### SKU B100-K200

- G-Tract A (Shroud, rump, etc.)
- G-Tract B (Vaginal Tube)
- G-Tract C (Uterine Horns)
- Rectum Bag & Clamp
- Blue Cervix or Red Cervix
- Tail
- Sim Poo
- Bungee & Carabiner Kit
- Blower & Blower Tube
- Pelvis Cushion
- Knurl Screws x 8

## 1b — INTERLOCKING BOVINE

### Mastitis Udder Kit

**SKU**  
**B100-K300**

- Vacuum Formed Udder
- Clear Teats x 4
- Mounting Rings x 4
- Teat Fasteners
- Knurl Screws x 4

### Milking Udder Kit

**SKU**  
**B100-K400**

- Vacuum Formed Udder
- Pink Teats x 4
- Mounting Rings x 4
- Teat Fasteners
- Knurl Screws x 4

### Metri- Treating Kit

**SKU**  
**B100-K500**

- G-Tract A (Shroud, rump, etc.)
- Uterine Horns with Stretchy V
- Rectum Bag & Clamp
- Bungee & Carabiner Kit
- Knurl Screws x 3

### Cupping Udder Kit

**SKU**  
**B100-K600**

- Cupping Udder on Mounting Board
- Knurl Screws x 4

# Bovine

## Calf Disbudding Kit

**SKU**  
**B300-K100**

**A model to teach cornual nerve block and horn disbudding techniques on calves.**

The Calf Disbudding Model incorporates an electronic cornual nerve block which allows you to practice feeling for the correct position of the nerve and injecting at the correct angle and depth. When the nerve is hit with the needle, a buzzer sounds and the LED light will change colour.

After the cornual nerve block has been mastered, you are able to practice disbudding. Using our Holsim Disbudding Tool and disposable horn buds, you can practice feeling

for the bud, cutting through the skin and flicking the bud out. The calf has a realistic feel with hard skull, soft neck and floppy reinforced ears. Our Calf Disbudding Model is battery powered so can be taken anywhere for class learning.

**KIT INCLUDES**

- 1 x Calf Disbudding Model
- 1 x Disbudding Tool
- 10 x Horn Buds

---

**Calf Disbudding Head Model ( B300-000 ) Available**

---





## Bovine

### Calf Disbudding Horn Buds Pack of 10

**SKU**  
B300-200

**Designed to be inserted into the Holsim Calf Disbudding Model.**

These disposable horn buds allow you to practise feeling for the horn bud, cutting through the skin and flicking the bud out.

A soft silicone skin encases the 3D printed horn and foam support ring. These are designed to be used with the Holsim Disbudding Tool. Standard pack includes 10 buds.



## Bovine

### Calf Disbudding Tool only

**SKU**  
B300-300

**Designed to be used with the Holsim Calf Disbudding Model and Horn Buds.**

The Holsim Disbudding Tool provides a realistic representation of a veterinarian's disbudding hot iron, used to disbud calves.

This tool allows you to practice the technique of burning through the skin and flicking the horn bud out, in a safe environment without the risk of burns or discomfort to the calf.



# Bovine

## Metri-Treating Model

**SKU**  
**B400-000**

**A model to teach students to detect cows with endometritis and treat uterine infections.**

Students are able to identify at-risk cows at calving time and treat any uterine infections with antibiotics. This teaches them to increase reproductive and economic performance.

Students should be able to recognise each body part of the cow through feel and diagnose endometritis with the detecting device.

The model has anatomically correct uterine and horns, vaginal tract, cervix, vulva and rectum of the cow. The Metri-Treating Model is a standalone table top set up for easy classroom learning and transportation.

The Metri-Treating Model packs down into separate components for easy compact travelling and set up is just as easy.

Silicone components and Holsim Sim Poo create a life-like tactile experience for the user.



# Bovine

## Milking Udder On Stand

**SKU**  
**B500-000**

**A model to teach milking and teat sealing techniques.**

Students can practise the technique of pulling and producing milk from the teats in a controlled environment with clear assessment from the tutor before advancing into the field. Students can also learn how treat mastitis in the teats.

The udder is made with moulded plastic, and each teat is removable for cleaning and maintenance. The udder can be filled with 'milk' for training.

There are two types of teats: pink for use with milking training, and clear, for teaching teat sealing techniques. The clear silicone allows students to see how far up into the teat the syringe is inserted and antibiotics are applied.



# Bovine

## Milking Teat

**SKU**  
**B500-K200**

**Designed to be used with the Holsim Milking Udder.**

These teats can be inserted into the milking udder for use of practicing milking techniques. Each teat can be cleaned after use. These teats come with mounting rings and fasteners for attaching to the milking udder.

**KIT INCLUDES**

4 x Pink Milking Teats

---

**Milking Teats only ( B500-200 ) Available**

---



## Bovine Mastitis Teat

**SKU**  
**B500-K300**

**Designed to be used with the Holsim Milking Udder.**

These teats can be inserted into the milking udder for use of practicing teat sealing for treatment of mastitis in cattle.

The clear silicone allows students to see how far up into the teat the syringe is inserted and antibiotics are applied.

Each teat can be cleaned after use. These teats come with mounting plates and fasteners for attaching to the milking udder.

**KIT INCLUDES**

4 x Clear Mastitis Teats

---

**Mastitis Teats only ( B500-300 ) Available**

---



# Bovine

## Udder Anatomy Model

**SKU**  
**B600-000**

**Cross-sectional anatomic display model of a cow udder.**

Holsim has created a silicone and foam section model of an udder for students to study and understand how the udder works in a controlled and evaluated environment.

Moulded foam and silicone create a realistic feeling udder, which is sectioned to reveal the mammary and teat structure.

The model replicates the mammary veins, arteries, lymphatic vessels, four quarters and milk passages as well as teat structure.

The udder can be taken off its stand for further investigation and learning in the classroom.



# Bovine

## Cow Tail with Blood Taking Kit

**SKU**  
**B700-K100**

**Model of a cow tail with artificial blood pump and life-like coccygeal vein for practicing blood taking technique.**

Vertebrae and vein positioning have been replicated from real bovine anatomy, allowing students to palpate for the correct position.

Blood is circulated through the vein by our Holsim Blood Pump to simulate the heart pumping. When inserting the needle, a realistic resistance is felt under the skin. Blood is circulated through the vein by our Holsim Blood Pump.

---

**Blood Sampling Tail & Insert ( B700-000 ) Available**

---



\* Note this does not come with clamp

When inserting the needle, a realistic resistance is felt under the skin. When pressure is applied, the vein bulges for students to see where to insert the needle.

The model is supplied with a sturdy metal stand which is to be clamped to a bench top, allowing the tail to hang down in a natural position.

### KIT INCLUDES

- 1 x Blood Sampling Tail with Coccygeal Vein Insert
- 1 x Blood Pump
- 1 x Jar of Artificial Blood
- 1 x Stand

*\*Replacement insert available*

AGRIHEALTH  
ULTRASOUND  
PRODUCTS



40 Day Pregnancy Palpation

|     |           |
|-----|-----------|
| SKU | B800-40PP |
|-----|-----------|

The Holsim 40-day ultrasound palpation uterine horn is filled with amniotic fluid and contains a 40-day fetus inside with small umbilical cord. The pregnant uterine horn has a slipped membrane inside, with visible creases and folds on the internal uterine wall. Placentomes are also visible and palpable on the internal uterine wall, enabling early-stage pregnancy diagnosis. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. The fetus has an air-filled brain inside the silicone body. This gives realistic imagery for pregnancy detection and diagnosis of the fetus. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.

70 Day Pregnancy Palpation

|     |           |
|-----|-----------|
| SKU | B800-70PP |
|-----|-----------|

The Holsim 70-day ultrasound palpation uterine horn is filled with amniotic fluid and contains a 70-day fetus inside with medium umbilical cord. The pregnant uterine horn has a slipped membrane inside, with visible creases and folds on the internal uterine wall. Placentomes are also visible and palpable on the internal uterine wall, enabling early-stage pregnancy diagnosis. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. The fetus has an air-filled brain and 3D printed rib and leg bones inside the silicone body. This gives realistic imagery for pregnancy detection and diagnosis of the fetus. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.

50 Day Pregnancy Palpation

|     |           |
|-----|-----------|
| SKU | B800-50PP |
|-----|-----------|

The Holsim 50-day ultrasound palpation uterine horn is filled with amniotic fluid and contains a 50-day fetus inside with small umbilical cord. The pregnant uterine horn has a slipped membrane inside, with visible creases and folds on the internal uterine wall. Placentomes are also visible and palpable on the internal uterine wall, enabling early-stage pregnancy diagnosis. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. The fetus has an air-filled brain and 3D printed rib bones inside the silicone body. This gives realistic imagery for pregnancy detection and diagnosis of the fetus. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.

80 Day Pregnancy Palpation

|     |           |
|-----|-----------|
| SKU | B800-80PP |
|-----|-----------|

The Holsim 80-day ultrasound palpation uterine horn is filled with amniotic fluid and contains a 80-day fetus inside with large umbilical cord. The pregnant uterine horn has a slipped membrane inside, with visible creases and folds on the internal uterine wall. Placentomes are also visible and palpable on the internal uterine wall, enabling early-stage pregnancy diagnosis. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. The fetus has an air-filled brain and 3D printed rib and leg bones inside the silicone body. This gives realistic imagery for pregnancy detection and diagnosis of the fetus. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.

60 Day Pregnancy Palpation

|     |           |
|-----|-----------|
| SKU | B800-60PP |
|-----|-----------|

The Holsim 60-day ultrasound palpation uterine horn is filled with amniotic fluid and contains a 60-day fetus inside with medium umbilical cord. The pregnant uterine horn has a slipped membrane inside, with visible creases and folds on the internal uterine wall. Placentomes are also visible and palpable on the internal uterine wall, enabling early-stage pregnancy diagnosis. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. The fetus has an air-filled brain and 3D printed rib bones inside the silicone body. This gives realistic imagery for pregnancy detection and diagnosis of the fetus. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.

90 Day Pregnancy Palpation

|     |           |
|-----|-----------|
| SKU | B800-90PP |
|-----|-----------|

The Holsim 90-day ultrasound palpation uterine horn is filled with amniotic fluid and contains a 90-day fetus inside with large umbilical cord. The pregnant uterine horn has a slipped membrane inside, with visible creases and folds on the internal uterine wall. Placentomes are also visible and palpable on the internal uterine wall, enabling early-stage pregnancy diagnosis. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. The fetus has an air-filled brain and 3D printed rib and leg bones inside the silicone body. This gives realistic imagery for pregnancy detection and diagnosis of the fetus. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.

50 Day Ultrasound Palpation

|     |           |
|-----|-----------|
| SKU | B800-50UP |
|-----|-----------|

The Holsim 50-day ultrasound palpation uterine horn is filled with amniotic fluid and contains a 50-day fetus inside with small umbilical cord. The pregnant uterine horn has a slipped membrane inside. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. The fetus has an air-filled brain and 3D printed rib bones inside the silicone body. This gives realistic imagery for pregnancy detection and diagnosis of the fetus. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.

80 Day Ultrasound Palpation

|     |           |
|-----|-----------|
| SKU | B800-80UP |
|-----|-----------|

The Holsim 80-day ultrasound palpation uterine horn is filled with amniotic fluid and contains a 80-day fetus inside with large umbilical cord. The pregnant uterine horn has a slipped membrane inside. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. The fetus has an air-filled brain and 3D printed rib and leg bones inside the silicone body. This gives realistic imagery for pregnancy detection and diagnosis of the fetus. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.



AGRIHEALTH  
ULTRASOUND  
PRODUCTS



60 Day Ultrasound Palpation

| SKU | B800-60UP |
|-----|-----------|
|-----|-----------|

The Holsim 60-day ultrasound palpation uterine horn is filled with amniotic fluid and contains a 60-day fetus inside with medium umbilical cord. The pregnant uterine horn has a slipped membrane inside. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. The fetus has an air-filled brain and 3D printed rib bones inside the silicone body. This gives realistic imagery for pregnancy detection and diagnosis of the fetus. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.

Empty/40 Day  
Pregnancy Palpation

| SKU | B800-40PE |
|-----|-----------|
|-----|-----------|

The Holsim Empty/40-day pregnancy palpation uterine horn is filled with amniotic fluid. Designed to be used with our Holsim AI Cow or Interlocking Bovine, the uterine horns lie on our Holsim Bean Bag inside the cow and can be easily splayed apart to lie flat for easy palpation. This gives realistic feel for pregnancy detection and diagnosis. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector with tabs, which plugs into our Holsim vagina tube. The uterine horns can be suspended by bungee cords through the tabs on the Uterine Y Body connector.

70 Day Ultrasound Palpation

| SKU | B800-70UP |
|-----|-----------|
|-----|-----------|

The Holsim 70-day ultrasound palpation uterine horn is filled with amniotic fluid and contains a 70-day fetus inside with medium umbilical cord. The pregnant uterine horn has a slipped membrane inside. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. The fetus has an air-filled brain and 3D printed rib and leg bones inside the silicone body. This gives realistic imagery for pregnancy detection and diagnosis of the fetus. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.

Uterine Y body W Tags

| SKU | B800-10PP |
|-----|-----------|
|-----|-----------|

The Holsim Uterine Y Body is designed to be used in conjunction with two of our Holsim Uterine Horns with Pregnancy Palpation. The Y body and two horns are cast together in one unit to prevent any amniotic fluid leaking. The Y Body then connects onto our Holsim Vagina Tube (G-Tract B.) The tags on the Y Body allow bungee cords to be connected, to be able to suspend the uterus inside our bovine models.

90 Day Ultrasound Palpation

| SKU | B800-90UP |
|-----|-----------|
|-----|-----------|

The Holsim 90-day ultrasound palpation uterine horn is filled with amniotic fluid and contains a 90-day fetus inside with large umbilical cord. The pregnant uterine horn has a slipped membrane inside. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. The fetus has an air-filled brain and 3D printed rib and leg bones inside the silicone body. This gives realistic imagery for pregnancy detection and diagnosis of the fetus. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.

40 Day Ultrasound Palpation

| SKU | B800-40UP |
|-----|-----------|
|-----|-----------|

The Holsim 40-day ultrasound palpation uterine horn is filled with amniotic fluid and contains a 40-day fetus inside with small umbilical cord. The pregnant uterine horn has a slipped membrane inside. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. The fetus has an air-filled brain inside the silicone body. This gives realistic imagery for pregnancy detection and diagnosis of the fetus. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.

Empty/40 Day  
Ultrasound Palpation

| SKU | B800-40UE |
|-----|-----------|
|-----|-----------|

The Holsim Empty/40-day ultrasound palpation uterine horn is filled with amniotic fluid. Designed to be used with our Holsim Ultrasound set up, the uterine horns have no air bubbles inside them, allowing for clear visual ultrasound imagery. This uterine horn is designed to be combined with another uterine horn and our Holsim Uterine Y Body connector without tabs, which plugs into our Holsim vagina tube.

Uterine Y Body without Tags

| SKU | B800-10UP |
|-----|-----------|
|-----|-----------|

The Holsim Uterine Y Body is designed to be used in conjunction with two of our Holsim Uterine Horns with Ultrasound Palpation. The Y body and two horns are cast together in one unit to prevent any amniotic fluid leaking. The Y Body then connects onto our Holsim Vagina Tube (G-Tract B.)



*\* The Rights of this product are owned by LIC. If you wish to purchase this must be done only through them.*

## 2a — LIC OWNED PRODUCTS

### Mk3 AI Cow

| SKU | B200-000 |
|-----|----------|
|-----|----------|

A transportable training model for teaching of artificial insemination techniques. The most realistic training device in the market replicates a cow's reproductive tract and pelvic cavity. It consists of a fibreglass body and silicone internals. Four windows allow internal views so the trainee's procedure can be guided and assessed. Heating simulates body temperature and inflation of the bowel mimics pneumorectum.

### AI Cow Arse ( G-Tract A)

| SKU | B200-200 |
|-----|----------|
|-----|----------|

Designed to attach to the Holsim AI Cow. This silicone rectum, vaginal entrance and rump are attached to a hard shroud which mounts onto the Holsim AI Cow backside using threaded knurl screws to secure it in place. The G-Tract B and C attach onto the vagina ring, allowing a smooth transition from vaginal entrance through the vagina tube with cervix to the uterus. The Holsim rectum bag attaches onto the rectum ring.

### Vagina Tube ( G-Tract B)

| SKU | B200-300 |
|-----|----------|
|-----|----------|

Designed to attach to the Holsim G-Tract A and G-Tract C for the AI Cow. This silicone tube is designed to replicate the vaginal canal of a cow. The Holsim Cervix can be inserted into the tube and moved along the vaginal tube to change the age of the cow. Both the Red (Medium) and Blue (Easy) Cervixes can be used in conjunction with this part of the G-Tract. We recommend using the Holsim Cervix Spoon to remove the Cervix after use.

### 60D Pregnant Uterus

| SKU | B800-90H60FE |
|-----|--------------|
|-----|--------------|

The Holsim 60-day pregnant uterine horn has a palpable silicone fetus and medium umbilical cord, while the other uterine horn is empty. The two uterine horns are filled with amniotic fluid, and can be easily splayed apart to lie flat on our Holsim Beanbag, for easy palpation. The pregnant uterine horn has a slipped membrane inside. The two uterine horns are connected with our Holsim Uterine Y Body connector, that plugs into our vagina tube. The uterine horns have tabs on the top of them to allow bungee cords to connect onto them, allowing them to be suspended inside the cow.

### Rectal Palpation Cyst

| SKU | B800-300 |
|-----|----------|
|-----|----------|

Designed to be used with our Pregnancy & Rectal Palpation Uterus, this interchangeable Ovarian Cyst is filled with soft foam, and attaches onto the end of the uterine horn. It enables students to palpate the ovaries to learn how to detect and diagnose various ovarian cysts during different cycle stages of the bovine.

### Rectal Palpation Ovary with CL

| SKU | B800-400 |
|-----|----------|
|-----|----------|

Designed to be used with our Pregnancy & Rectal Palpation Uterus, this interchangeable Ovary with CL has a large concave end, and attaches onto the end of the uterine horn. It enables students to palpate the ovaries to learn how to detect and diagnose ovaries with CL during different cycle stages of the bovine.

### Cow Cervix (Blue - Easy)

| SKU | B200-500 |
|-----|----------|
|-----|----------|

Designed to be used in the Holsim Vagina Tube. This silicone cervix is used for artificial insemination training with the Holsim AI cows. The cervix displays a large cervical opening. This cervix is rated as an easy level of insemination difficulty. Rings inside this cervix represent the annular folds. These allow students to learn how to manipulate these rectally when an insemination rod is passed through to the uterus.

### Cow Cervix ( Red - Medium )

| SKU | B200-600 |
|-----|----------|
|-----|----------|

Designed to be used in the Holsim Vagina Tube. This silicone cervix is used for artificial insemination training with the Holsim AI cows. The cervix displays a smaller cervical opening than the blue cervix. This cervix is of medium difficulty to inseminate. Rings inside this cervix represent the annular folds. These allow students to learn how to manipulate these rectally when an insemination rod is passed through to the uterus.

### Uterine horns ( G-Tract C )

| SKU | B200-400 |
|-----|----------|
|-----|----------|

Designed to attach onto the Holsim G-Tract-B. These are silicone uterine horns, which are the third part of the Holsim G-Tract. Designed to be used in the Holsim AI cow, students can learn to palpate the uterus and ensure sperm is applied evenly into the uterine horns for effective and successful artificial insemination. Ovaries can be palpated at the end of each uterine horn.

### Rectal Palpation Ovary — Follicle

| SKU | B800-500 |
|-----|----------|
|-----|----------|

Designed to be used with our Pregnancy & Rectal Palpation Uterus, this interchangeable Ovary with follicle has a soft hollow end, and attaches onto the end of the uterine horn. It enables students to palpate the ovaries to learn how to detect and diagnose ovaries with follicles during different cycle stages of the bovine.

### Rectal Palpation Ovary — Polyps

| SKU | B800-600 |
|-----|----------|
|-----|----------|

Designed to be used with our Pregnancy & Rectal Palpation Uterus, this interchangeable Ovary with small follicles/polyps attaches onto the end of the uterine horn. It enables students to palpate the ovaries to learn how to detect and diagnose small follicles or polyps on the ovary during different cycle stages of the bovine.



2a —

## LIC OWNED PRODUCTS

### BOVINE

#### Rectal Palpation & Pregnancy Kit

##### Teaching aid for manual pregnancy check and ovarian diagnostics.

Students are now able to practice the most important steps of manual pregnancy checks and ovarian diagnostics.


Placing the uterus on an intestine cushion/bean bag, allows for realistic feeling and students are able to palpate the entire tract, including cervix, uterine horns and different variations of ovaries."

*\* The Rights of this product are owned by LIC.  
If you wish to purchase this must be done  
only through them.*

important signs of an early pregnancy (enlarged, fluid filled horn and membrane slip.) Additionally, the trainer can set up different cycle stages by exchanging ovaries which feature different functional bodies (small follicles, large follicle, Corpus luteum, cyst.)

This model represents the uterus with one horn in a pregnancy stage of approx. 42 days, and the other at 60 days. The pregnant horn is enlarged with a foetus contained inside. The pregnant horn is filled with amniotic fluid. Two membranes enable training of manual membrane slipping for pregnancy diagnostics.

\* Spare parts available

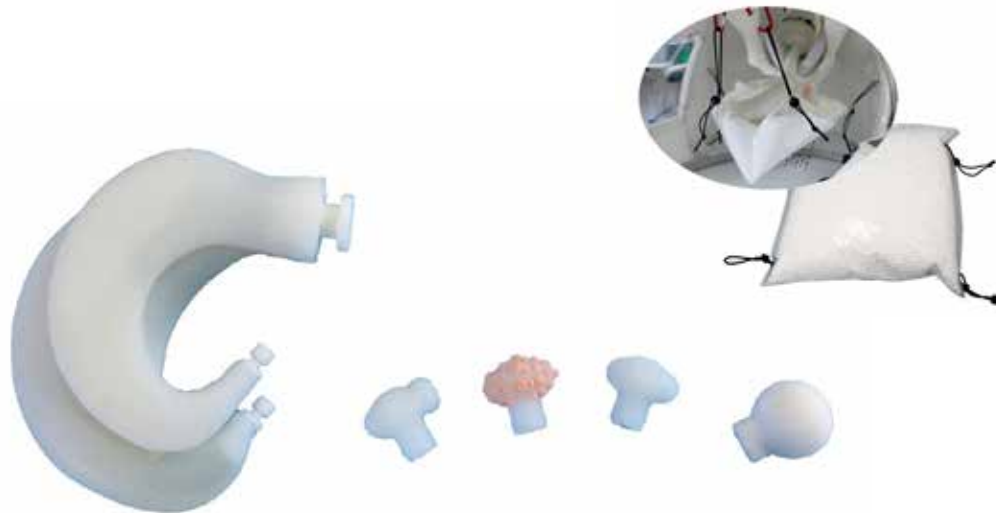

*There is lower stress on the student and lower emotional and psychological effects when performing operations.*

*Training on models decreases the risk of future vets becoming desensitised to animals being in pain.*



03 —

# Canine



# Canine

## Dog Catheter Model ( Female )

**SKU**  
**C100-000**

**Section of back end of dog - approximately Labrador sized. Vagina, clitoral fossa, urethra and bladder in life-like cast silicone.**

Mesh reinforced vulva for durability, a section of pelvis bone and soft foam surrounds the silicone g-tract for a life-like feel.

The bladder is filled with water from a tube at the back, which clips out of the way when in use. Once a catheter is successfully inserted, this water can be drained from the bladder like it would be from a real dog.

The model can be used upright (as shown) or on its side, with a bench mount for stability in both orientations.



# Canine

## Brachycephalic Dog Model

**SKU**  
**C200-000**

**A model to teach intubation techniques on Brachycephalic dogs.**

The Brachycephalic Intubation Model helps you navigate the mouth and respiratory region of Brachycephalic (short muzzle) dogs to help improve your ability to perform an intubation operation in real-life situations.

This anatomically correct training aid features a functioning jaw, a stretchy tongue, and an accurate throat system (epiglottis, trachea, soft palate, etc.). The throat system is visible at the back so you can check your accuracy in performing the procedure.

The intubation model also comes with 4 mounting brackets to mount to a stand or other device of your choice.



# Canine

## Canine Intubation Model

**SKU**  
**C300-000**

**A model to teach intubation techniques on medium to large sized, long muzzle dogs.**

The Canine Intubation Model helps you navigate the mouth and respiratory region of a dog to help improve your ability to perform an intubation operation in real-life situations.

This anatomically correct training aid features a functioning jaw, a stretchy tongue, and an accurate throat system (epiglottis, trachea, soft palate, etc.) The throat system is visible at the back so you can check your accuracy in performing the procedure.

The intubation model also comes with 4 mounting brackets to mount to a stand or other device of your choice.



# Canine

## Oesophageal Dog Model

**SKU**  
**C400-000**

**A training model for classroom teaching of feeding tube insertion into dog oesophaguses.**

The Canine Oesophageal Feeding Tube model helps you learn where to make an incision on a dog's neck and how to navigate tools through the mouth and respiratory region to improve your ability to place a feeding tube through a dog's neck.

This training aid features a functioning jaw, accurate mouth and throat features (tongue, epiglottis, soft palate, etc.), and a surgical block with a jugular vein.

The surgical jugular block is a replaceable part that can be cut and sutured multiple times. It also allows you to palpate where the jugular vein is located to improve your operational skills and connect the blood pump to the jugular vein.

There is a stomach opening at the rear of the model to show you if you have correctly placed your feeding tube. The model comes with connection points on the back to attach a mounting bracket if needed. Contact us for any spare parts or mounting brackets.





## Canine

### Canine Oesophageal Jugular Block

**SKU**  
**C400-200**

**Replacement jugular block for the Holsim Canine  
Oesophageal Feeding Tube Model.**

The surgical jugular block can be cut and sutured multiple times. It also allows you to palpate the jugular vein, and learn where it is located, to improve your operational skills.

The jugular block sits into the dog model snugly, intersecting the trachea and oesophagus. The block is secured in place with velcro.



# Canine

## Ear Examination Model ( Otoscope )

**SKU**  
C500-000

Ear problems in dogs are a frequent occurrence in veterinary clinics. Most cases can be resolved with a thorough examination and clean of the ear. Ear infections are often caused by allergy-related issues, foreign bodies or the anatomy of the ear.

The ear model includes an anatomically correct silicone pinna, representing that of a Labrador, and an ear canal (horizontal and vertical). The pinna can be pulled gently up and back from the head of the dog, to open up the canal tract.

Students can learn how and where to palpate the ear canal correctly to dislodge any debris trapped in the ear. The palpable area is made from softer silicone than the pinna, allowing it to be pushed and rubbed in a circular motion to dislodge any debris.

Due to the model being made from silicone, the ear canal can be filled with cleaning fluid or saline solution.

The pinna is detailed with cartilage folds, to allow students to practice cleaning the outer area of the pinna.

The ear canal can also be filled with substances to represent wax or debris in the ear, giving a more realistic representation of the ear. The ear can be easily cleaned and flushed out with soapy water after use.

Interchangeable coloured ear drums allow the tutor to challenge students on what to look for when using the otoscope down the ear canal. These ear drums come as a set of four "ear drum plugs." One is placed at the end of the horizontal canal during the procedure.

The ear models are mounted to an interlocking plywood dog head. The normal ear block has a pale pink inner ear and canal with a wide opening. The infected ear block has a red swollen entrance to the ear and a narrower and stiffer ear canal.



# Canine

## Jugular Blood Sampling Kit

**SKU**  
**C600-K100**

**Canine head and neck model with jugular vein and blood pump.**

Anatomically correct head and neck of a medium sized dog, with replaceable jugular vein patch on one side.

Vein patch has skin, fat and muscle layers to give a more realistic feel when inserting needle. Blood is circulated through the vein by our Holsim Blood Pump to simulate the heart pumping.

When inserting the needle, a realistic resistance is felt under the skin. When pressure is applied, the vein bulges for students to see where to insert the needle.

**KIT INCLUDES**

- 1 x Canine Jugular Blood Sampling Model
- 1 x Blood Pump
- 1 x Jar of Artificial Blood

---

**Jugular Blood Sampling Model (C600-000) Available**

---



# Canine

## Canine Cardio Audio Box

**SKU**  
**C700-000**

**A training model to identify heart murmurs in each heart valve in dogs.**

Veterinary students are able to use their stethoscope to determine which of the four heart valves have murmurs. The audio box enables the tutor to associate murmurs with individual heart valves for the student to listen to and identify. The Canine Cardio Audio Box has been designed to only be heard with the stethoscope during student training when inserted into a Holsim Dog.

---

**With Lung Auscultation ( C700-200 ) Available**

---

The audio box anatomically mirrors the heart of a dog with individual murmurs to be associated with each heart valve. Students can find the heart placement by counting the ribs of the dog.

Each valve can be turned on or off to include a "whoosing" sound to represent the murmur. The audio box represents the four heart valves: mitral, aortic, pulmonary and tricuspid. The box is able to undergo heart compressions when inserted into a CPR Dog. An additional module to allow lung auscultation training with a variety of lung sounds can also be added.

*Note: A CPR Dog must be provided in conjunction with these products. We can provide this with extra cost added, for more information get in touch with us via email.*



# Canine

## IV Dog Leg Kit Standard

**SKU**  
**C800-K100**

**A simple, transportable training tool for fluid administration. This is a kit that includes the model, blood pump, IV insert and jar of blood.**

Veterinary students are able to practise finding the cephalic vein in a controlled environment for intravenous therapy or blood sampling with clear assessment and feedback from the tutor.

The dog skin insert is firstly connected to the tubes in the dog leg. The blood is poured into the pump and then switched on. The blood will circulate through the vein tubing, where the student can then attempt to identify the cephalic

---

**IV Leg Model Only ( C800-000 ) Available**

---

vein and withdraw blood or administer the required fluids or drugs. The Holsim IV Dog Leg allows students to understand the bone and vein structure of a dog leg. The shaved skin can be removed to show the layer of skin, fat and muscle tissue. The blood is pumped through the main cephalic vein into the skin section where students, if accurate, will be able to withdraw blood, inject drugs or practice setting up a drip.

### **KIT INCLUDES**

- 1 x IV Dog Leg Model with IV Insert
- 1 x Blood Pump
- 1 x Jar of Artificial Blood



# Canine

## IV Dog Leg Insert – Standard

**SKU**  
**C800-200**

**Designed to be used with the Holsim IV Dog Leg and Canine Jugular Blood Sampling Model.**

The dog skin insert is first connected to the tubes in the dog leg or neck. The blood will circulate through the vein tubing, where the student can then attempt to identify the vein and withdraw blood or administer the required fluids or drugs.

The insert has a realistic resistance under the skin when the needle is inserted. The vein also bulges slightly when pressure is applied. The insert also has coloured skin, fat and muscle layers.



# Canine

## Hind Limb IV Bandaging Leg Kit

**SKU**  
**C800-K1100H**

**An IV and Bandaging training model all in one. This is the whole kit including range of products such as Blood pump, Artificial blood and IV Dog Limb Model with hind limb insert.**

Veterinary students are able to practice finding the cephalic vein in a controlled environment for intravenous therapy or blood sampling with clear assessment and feedback from the tutor.

---

**IV & Bandaging Hind Limb ( C800-1000H ) Available**

---

Blood circulates from the Holsim Blood Pump through the vein tubing, where the student can then attempt to identify the cephalic vein and withdraw blood or administer the required fluids or drugs.

The leg can also be used for teaching bandaging methods, especially as the toes of the dog leg can be spread to wrap the bandage between.

### **KIT INCLUDES**

- 1 x IV Dog Hind Limb Model with IV Insert
- 1 x Blood Pump
- 1 x Jar of Artificial Blood



# Canine

## IV Insert - Bandaging Hind Limb

**SKU**  
**C800-1200H**

**Designed to be used with the Holsim  
Hind Limb Bandaging & IV Model.**

The dog skin insert is first connected to the tubes in the dog leg. The blood will circulate through the vein tubing, where the student can then attempt to identify the vein and withdraw blood or administer the required fluids/drugs.

The insert has a realistic resistance under the skin when the needle is inserted. The vein bulges slightly when pressure is applied. The insert also has coloured skin, and muscle layers.





# Canine

## Forelimb Canine Bandaging & IV Kit

**SKU**  
**C800-K1100F**

### **An IV and Bandaging training model all in one.**

Veterinary students can practise finding the cephalic vein in a controlled environment for intravenous therapy or blood sampling with clear assessment and feedback from the tutor.

Blood circulates from the Holsim Blood Pump through the vein tubing, where the student can then attempt to identify the cephalic vein and withdraw blood or administer the required fluids/drugs.

The leg can also be used for teaching bandaging methods, especially as the toes of the dog leg can be spread to wrap the bandage.

### **KIT INCLUDES**

- 1 x IV & Bandaging Forelimb Model with IV Insert
- 1 x Blood Pump
- 1 x Jar of Artificial Blood

---

**IV & Bandaging Forelimb ( C800-1000F ) Available**

---



# Canine

## IV Insert - Forelimb Bandaging Model

**SKU**  
**C800-1200F**

**Designed to be used with the Holsim Forelimb Bandaging  
& IV Model**

The dog skin insert is first connected to the tubes in the dog leg. The blood will circulate through the vein tubing, where the student can then attempt to identify the cephalic vein and withdraw blood or administer the required fluids/drugs.

The insert has a realistic resistance under the skin when the needle is inserted. The vein bulges slightly when pressure is applied. The insert also has coloured skin, and muscle layers.



# Canine

## Articulated Dog Bones

**SKU**  
**C900-200**

**Anatomically precise 3D printed dog bones.**

3D printed bones accurately show important osteological features of the hind limb, from the hip down to the digits, enabling students to learn and identify the name, position and function of the features.

The bones come pre-drilled in exactly the right places (for threading together) to allow easy demonstration of how the joints fit together and the major ligaments of joints. Bones can be drawn on with water-based pens and then washed clean in soapy water.

Holsim Dog Bones were 3D scanned from bones prepared from a single real dog, and then 3D printed from renewable bioplastic. Because the original scan data for the bones came from a single dog, all of the bones are in proportion and fit together snugly.

This product can be purchased in both pre assembled as well as disassembled.

**Disassembled ( C900-000 ) or Assembled Available**



# Canine

## Dog Bones

**SKU**  
**C1000-000**

**Anatomically precise 3D printed dog bones enable students to name, position and function each bone and its features.**

Holsim Dog Bones were 3D scanned from bones prepared from a single real dog, and then 3D printed from bioplastic. Because the original scan data for the bones came from a single dog, all of the bones are in proportion and fit together snugly.

The fore and hind limbs have all major bones, all carpals and tarsals and a representative digit (III). The set includes vertebrae from each region, a sacrum and two ribs.

Limb bones in the Holsim Dog Bones set come pre-drilled for threading together, demonstrating how the joints fit together with the major ligaments.



# Canine

## Mesopcephalic Dental Dog

**SKU**  
**C1100-000**

**Holsim's Mesocephalic Dental Dog is an anatomically correct training model for teaching dental and intubation techniques on medium-sized dogs.**

Our canine dental training model features a full set of replaceable resin teeth, with soft silicone gums. The teeth include all 42 teeth, with incisors, canines, premolars, and molars. Each section can be swapped in or out easily with a screwdriver.

The dental model features an open jaw. The dental model features a stretchy silicone tongue, soft palate, epiglottis, and entrances to the trachea and oesophagus for teaching intubation techniques.



# Canine

## Mesocephalic Dental Dog - Teeth Set (Standard)

**SKU**  
**C1100-200**

**Holsim Mesocephalic Standard Teeth Set is designed to be used with our Mesocephalic Dental Dog Training Model.**

The teeth are 3D printed from resin and have a soft pink silicone gum attached. The set includes all 42 adult teeth, with incisors, canines, premolars, and molars for both the upper and lower jaw.

The teeth are divided into sections to be able to be swapped in and out of our dental dog head easily with a screwdriver.

Ask us about different types of replacement teeth sets! We will be having a tooth extraction set, and plaque-covered teeth set in the future. Stay tuned!



*Societal attitudes toward educational animal use, particularly with invasive procedures, have changed in recent decades. Regulations regarding animal welfare have also changed, making the traditional uses of animals in the teaching of clinical procedures less acceptable from both societal and regulatory perspectives.*



04 —

# Equine



# Equine

## Distal Limb Nerve Block

**SKU**  
**E100-000**

**An easy-to-use model that trains students effectively in distal nerve block procedures.**

An anatomically correct accurate representation of a horse's leg below the carpus.

*Students can practice the following nerve blocks:*

- Low Palmar digital
- Mid-Pastern Palmar digital
- Low Palmar
- Palmar

Audio feedback is sounded upon the correct placement of nerve block needles. The model has life-like nerve palpation and needle resistance using our electro-silicone technology.

---

**Distal Limb Insert (E100-200) Available**

---





# Equine

## Nasogastric Intubation Horse

**SKU**  
**E200-000**

**Teaching aid for the inserting and placement of nasogastric tubes through the oesophagus.**

Holsim's Equine Nasogastric Intubation Model is an anatomically correct teaching aid for veterinary students to learn how to insert nasogastric tubes through the oesophagus to relieve colic issues, administer drugs, or provide nutrition to the horse when not eating. This teaching aid can also be used to teach tube entry into the trachea for respiratory issues. Nasogastric intubation is an essential and possibly life-saving procedure performed routinely in cases of equine colic to decompress the stomach and provide therapy.

The model has a soft silicone nose, with both true and false nostrils. The nasal cavities lead up through a 3D printed resin skull, representing the hard palette in the nasal cavity. As the tube reaches the pharynx, a soft resistance can be felt as it connects with the larynx/oesophagal opening. The head of the horse can either be pushed downwards for the tube to pass into the oesophagus or can be raised for the tube to enter the trachea. The tube can be turned and adjusted both

ventrally and medially to ensure correct positioning before entering either passage. The tube can visually be seen travelling down the left-hand side of the horse's neck, as it passes down through the oesophagus. Students can also physically feel the tube passing down underneath the horse's skin on the neck

The trachea has a slight ribbing to it, to represent the tracheal rings/cartilages. These can be felt when the tube is passed down the trachea. If the tube is in the trachea, shaking it will result in a hollow rattle.

The Holsim Equine Nasogastric Model can be used in the classroom for ease of teaching. The base of the model has a drainage hole where a stomach can be connected underneath, filled with sour feed or stomach contents to teach students syphoning of the stomach. The neck skin on the horse can be removed for washing.



# Equine

## Rectal Palpation and Ultrasound Training Model

**SKU**  
**E300-000**

The Holsim Equine Rectal Palpation and Ultrasound Training Model is a state-of-the-art teaching tool designed for veterinary students and professionals to practice essential equine reproductive examination techniques. This life-size model replicates the back half of a horse with accurate anatomical structures, allowing for hands-on training in rectal palpation, ultrasound imaging, and twin pregnancy management.

### Key Features

**Life-Size Construction:** Built with interlocking plywood and mounted on a sturdy base, the model provides a realistic and durable training experience.

**Anatomically Accurate Design:** Includes a hard pelvis, a fibreglass rump with a horsehair tail, and a realistic reproductive system made from high-quality silicone.

**Detailed Reproductive Anatomy:** Features a uterus, cervix, ovaries, and vaginal canal crafted from silicone for lifelike palpation and examination.



**Twin Pregnancy Simulation:** The uterine horns have bifurcations. One uterine horn contains a single 13-day pregnancy, while the other uterine horn contains twin 13-day pregnancies. Additionally, the main uterus body includes a 25-day pregnancy, allowing students to detect and diagnose different stages and complexities of pregnancy.

**Realistic Ultrasound Imaging:** Produces clear and accurate ultrasound visuals, helping students refine their imaging techniques and differentiate pregnancies from cysts.

**Palpation Practice:** The G-Tract assembly rests on a soft beanbag, allowing students to practice nuanced palpation techniques, including the identification and management of twin pregnancies.

**Rectal Training:** Entrance through a silicone rectum into a rectal bag enables students to practice proper palpation techniques safely.



EQUINE

**Instructor-Friendly Design:** Allows tutors to easily demonstrate procedures, assess student performance, and provide targeted feedback.

**Easy Setup and Portability:** Designed for quick assembly and disassembly, making it convenient to transport in a car, ute, or van for mobile training sessions.

### Addressing Industry Needs:

Veterinary professionals and students require extensive training to accurately detect and manage equine pregnancies, particularly twin pregnancies, which pose significant risks to both mares and foals. This model provides an invaluable solution for:

- Training vets to differentiate pregnancies from cysts and other uterine structures.
- Improving ultrasound technique for clearer, more accurate imaging.
- Enhancing manual palpation skills to detect and manage pregnancies.
- Practicing twin pregnancy separation and 'crushing' procedures in a controlled, low-risk environment.



05 —

# Feline

*For many students, participation in harmful live animal use, generates powerful emotional experiences and high levels of stress. These have considerable potential to adversely affect cognitive processes such as learning.*

# Feline

## Male Urinary Catheter

**SKU**  
**F100-000**

**Male feline model that allows practice of correct urinary catheter placement.**

This model is an anatomically correct male feline that allows practice of correct urinary catheter placement. The model can be clamped to a desktop or table and has non-slip rubber feet to prevent it from moving during the procedure. An internal casing can be removed from the cat, which holds the bladder, pelvis bone and penis with a soft foam surrounding.

A replaceable suture block sits in the backside of the cat, with a mesh reinforcement in the silicone, allowing the catheter to be stitched to the skin. The bladder can be filled with water from the fill tube on the back of the stand, for real life feedback and bladder drainage practice. The prepuce can be rolled backwards to expose the penis for catheter insertion.

We have a range of colour options for the different cats such as grey&white, ginger, black, and black&white.



# Feline

## Male Feline Urinary Catheter Suture Block

**SKU**  
**F100-200**

**Replacement suture block designed to be used with the Holsim Male Feline Urinary Catheter Model.**

This is the suture block for the Male Feline Urinary Catheter Model. This suture block sits in the backside of the cat and snaps into place magnetically.

The suture block has mesh reinforcement in the silicone, allowing the catheter to be stitched to the skin. The block includes anus and prepuce.





06 —

# Porcine

*Approximately 90% of published educational evaluations have shown that students being taught with ethical teaching methods achieved the same or better learning outcomes than using live animals.*

# Porcine

## Interlocking Plywood Pig


**SKU**  
**P200-000**

**A flat-pack model allows for a variety of veterinary models to be used in conjunction with it.**

This life-size sow is a flat-pack, interlocking model which allows feature models to slot into compartments on the body and head. This allows you to upgrade the pig with new procedures as they're developed. Made from high-grade 12mm ply with a coating/finish of your choice, this pig is robust and durable for student training in the classroom.

The head can raise and lower 30 degrees to the horizontal, for use when blood sampling from the jugular and restraining/snaring the pig's snout. When the head lowers, the spring-loaded jaw will open. When the head is raised, the jaw will shut closed.

Durable 12mm high-grade plywood frame. Flat-pack for ease of transportation/shipping. Articulating neck and spring-loaded jaw for raising/lowering head to +/- 30degrees. Able to have Holsim feature models added/removed at ease. Removable top cage for ease of access to G-tract set up. Life-size anatomically correct sow.

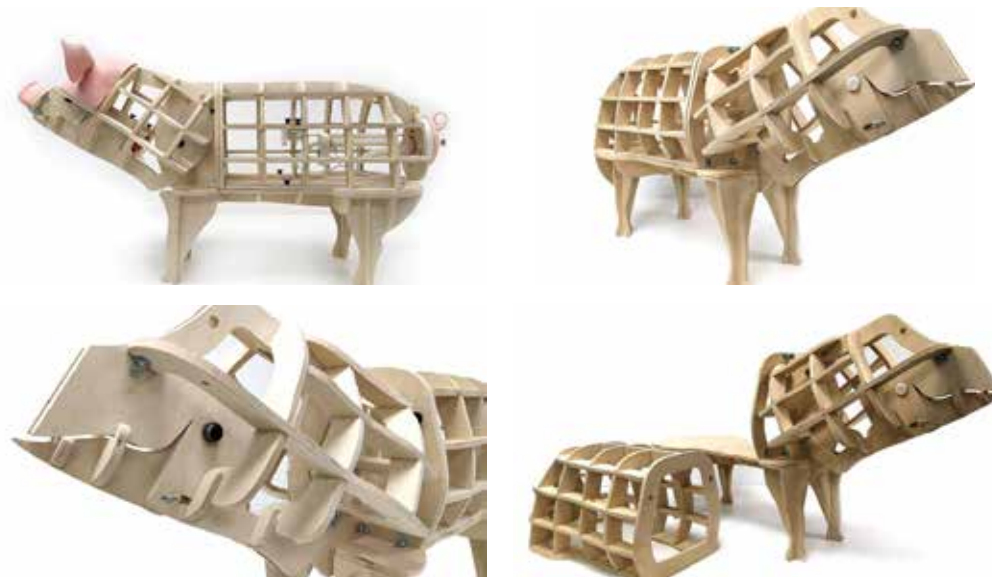



# Porcine

## Intramuscular Block With Blood Taking Kit

**SKU**  
**P300-K100**

**A training model designed to be used with the Holsim Interlocking Porcine for Intramuscular Injection training and Auricular Vein Sampling.**

This anatomically correct intramuscular injection block has conductive silicone patches embedded under the skin behind the ear. When the correct injection point has been achieved, a buzzer sounds, and the LED light turns green.



There is also an option to have the buzzer off while the student is training, when they believe they're at the right point of injection, the tutor may switch the buzzer on to see if the student is correct.

This anatomically correct pig ear model has a palpable silicone lateral auricular vein running through the ear which connects to our Holsim Blood Pump. The ear can be folded to expose the vein for needle insertion.

### **KIT INCLUDES**

- 1 x Intramuscular Injection with Blood Sampling Ear
- 1 x Blood Pump
- 1 x Jar of Artificial Blood
- 1 x Knurl Screw



# Porcine

## Pig Ear

**SKU**  
**P400-000**

**A training model is designed for lateral auricular vein blood sampling kit.**

This anatomically correct pig ear model has a palpable silicone lateral auricular vein running through the ear which connects to our Holsim blood Pump. The ear can be folded to expose the vein for needle insertion. This ear is stand-alone, however for a full experience purchase the Intramuscular Block & Blood Sampling Ear combination.



# Porcine

## Ear Blood Taking Kit

**SKU**  
**P400-K100**

**A training model for lateral auricular vein blood sampling.**

This anatomically correct pig ear model has a palpable silicone lateral auricular vein running through the ear which connects to our Holsim Blood Pump. The ear can be folded to expose the vein for needle insertion.

This ear is stand-alone, however, for a full experience, purchase the Intramuscular Block & Blood Sampling Ear Combination.

### KIT INCLUDES

- 1 x Blood Sampling Ear
- 1 x Blood Pump
- 1 x Jar of Artificial Blood



# Porcine

## Snout Nose Ringing/ Snaring Model

**SKU**  
**P500-000**

**A soft silicone snout with a tough septum for nose ringing and snaring techniques.**

Tough cartilage is embedded in the septum and tip of the snout for various ringing application training techniques

The snout is made from soft silicone. This can be bought stand-alone or can be bought for mounting onto our Interlocking Porcine Model for a more life-like experience and can be used with the articulating neck of the pig and snaring techniques.



PORCINE

# Porcine

## Jugular Blood Taking Kit

**SKU**  
**P600-K100**

**A training model for teaching jugular blood sampling, for use with the Interlocking Porcine.**

This jugular block has a definite visual curve, following the hollow of the neck, to indicate the injection point for blood sampling from the jugular vein.

Soft foam is embedded under silicone skin for a tender feel during palpation. Our Holsim Blood Pump connects to the jugular vein to allow for a pulsing blood flow.

This jugular block is best used with our Interlocking Porcine Model as the head can be raised to 30 degrees from the horizontal to correctly expose the jugular vein for blood sampling.

### KIT INCLUDES

- 1 x Jugular Blood Sampling Model
- 1 x Blood Pump
- 1 x Jar of Artificial Blood
- 2 x Knurl Screw

---

**Jugular Blood Sampling Model (P600-000) Available**

---



# Porcine

## Interlocking AI Model

**SKU**  
**P700-K100**

**A training model for teaching artificial insemination, for use with the Interlocking Porcine.**

Anatomically correct silicone internal tract with the cervix, vaginal tube and bladder entrance. When the correct place is hit with the insemination rod, an "oink" is sounded. Practice makes perfect because if you get it wrong, it hurts! A badly inseminated sow won't let anyone inseminate her again, forcing your prized sow to become bacon. This is best when mounted to our Interlocking Porcine for a more realistic experience.

**KIT INCLUDES**

1 x Artificial Insemination Model  
1 x Keyhole Bracket  
1 x Shroud  
2 x Knurl Screw

# Accessories

*Regarding practical AI training, the aspect of animal welfare in modern animal breeding is currently moving more and more into the focus. A study was carried out at the German Institute for Reproduction of Farm Animals IFN in cooperation with the University of Veterinary Medicine in Hanover to assess stress levels of cattle during artificial insemination. It was found that participants who had carried out preliminary exercises on our Holsim AI Cow model, induced lower cortisol concentrations in the cows than participants who were only theoretically instructed beforehand. Therefore, the research group recommends instructing the participants on an animal model before the practical exercise.*

# Accessories

## Suture Pad

**SKU**  
**Z100-000**

### A training pad for practising surgical stitching.

Holsim Suture Training Pads allow students to maximise surgical stitching techniques. Tutors can teach their students needle placement, stitching patterns, knot-tying, superficial and surgical closures.

Suture Training Pads can come in a range of wounds, cut patterns and skin layer depths. Students will learn how to select the correct needle and thread type to use for the wound, needle insertion positioning, stitch spacing, burying the final knot and removing the stitches once the wound has healed.

Holsim Suture Pads will enhance student learning and allow them to develop their skills, prepare the wound for sutures, understand the properties of suture equipment, and gain a full understanding of skin anatomy.

The suture pad replicates the skin and tissue layers; epidermis, subcutaneous and deep fascia for a realistic learning experience. Skin layers are reinforced with mesh to give a lifelike tension and strength. Suture pads provide a realistic resistance for the needle to push through during stitching.



*\*Note these images are the older version has been slightly upgraded*

## Accessories

### Suture Pad Base

**SKU**  
**Z100-200**

**A non-slip base for mounting the Holsim Suture Pad.**

The Holsim Suture Pad non-slip base fits our standard suture pads, raises the pad off the table, and provides a firm non-slip surface for practising suturing.

The non-slip base comes in white, with other colours are available for larger orders.

Contact us for bulk order pricing.



*\*Note these images are the older version has been slightly upgraded*



# Accessories

## Suture Kit


### SKU

Z100-K100

#### A training kit with suture pad and non-slip base for practising surgical stitching.

Holsim Suture Training Pads allow students to maximise surgical stitching techniques. Tutors can teach their students needle placement, stitching patterns, knot-tying, superficial and surgical closures.

Suture Training Pads can come in a range of wounds, cut patterns, and skin layer depths. Students will learn how to select the correct needle and thread type to use for the wound, needle insertion positioning, stitch spacing, burying the final knot, and removal of the stitches once the wound has healed.

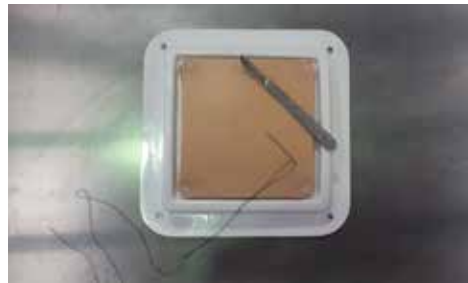

*\*Note these images are the older version has been slightly upgraded*

Holsim Suture Pads will enhance student learning and allow them to develop their skills, prepare the wound for sutures, understand the properties of suture equipment, and gain a full understanding of skin anatomy. Replicates the skin and tissue layers; epidermis, subcutaneous and deep fascia for a realistic learning experience. Skin layers are reinforced with mesh to give a lifelike tension and strength. Suture pads provide a realistic resistance for the needle to push through during stitching.

A non-slip base for mounting the Holsim Suture Pad. The Holsim Suture Pad non-slip base fits our standard suture pads. It raises the pad off the table and provides a firm non-slip surface for suture practice.

#### KIT INCLUDES

- 1 x Suture Pad
- 1 x Non-Slip Base
- 4 x Push Fasteners

## Accessories

### Amniotic Fluid

**SKU**  
**Z200-000**

**Holsim Amniotic Fluid is designed to be used with the Holsim Interlocking Bovine and the Holsim Birthing Calf.**

When practising the birthing process, the Birthing Calf should be covered in Amniotic Fluid while being placed in the uterus bag. This represents the amniotic sac as realistically as possible and allows the calf to slide out of the birth canal naturally and smoothly.



# Accessories

## Sim Poo

**SKU**  
**Z300-000**

**Holsim Sim Poo is designed to be used with the Holsim Interlocking Bovine.**

Sim Poo provides the cow with a realistic representation of excrement inside the rectum, giving the students an accurate tactile experience during AI training.

The Artificial Insemination Kit includes a barreling pump to mimic the action of pneumonectomy in cows, and the Sim Poo is used to teach students how to relieve rectum pressure, to avoid being sprayed with Sim Poo.



# Accessories

## Blood Pump

**SKU**  
**Z400-000**

**Designed to be used with the Holsim Blood Sampling Models.**

This is a small battery-powered, electronic pump, that replicates the rhythm of a heart pumping.

The pump will circulate Holsim Artificial Blood through the vein tubing and into the veterinary model e.g. IV or Jugular patches. The blood is pumped through vein tubing into the blood sampling model where students, if accurate, will be able to withdraw blood, inject drugs or practice setting up a drip. Available continuous flow for tail vein



*\*Note these images are the older version  
has been slightly upgraded*

# Accessories

## Artificial Blood

**SKU**  
**Z400-200**

**Designed to be used with the Holsim Blood Pump and Holsim blood sampling veterinary models.**

This artificial blood simulates the colour and consistency of real blood. Students can use this to practice blood sampling techniques accurately. The jar contains 400ml of artificial blood.



## Accessories

### Blower & Blower Tube

**SKU**  
**Z500-000**

**Designed to be used with the Holsim Inflatable Bag for birthing simulation.**

This is an electronic pump for filling the inflatable bag to support the calf in the uterus bag, providing a cushioned surrounding. The blower can plug into the wall for power for fast inflation.



## Accessories

### Bungee & Carabiner Kit

**SKU**  
**Z700-000**

**Designed for a variety of functions across the Holsim product range.**

In the Bovine Artificial Insemination Kit, the bungee cords are used to adjust the position of the uterus and cervix to simulate different ages of cows.

In the Bovine Birthing Kit, they are used to suspend the uterus bag to make it easier or harder to birth the calf. It is also used to allow the simulation of a prolapsed uterus.

They can be used to adjust the tension of the cervix to simulate degrees of difficulty.



## Accessories

### Holsim Wooden Mallet

**SKU**  
**Z900-000**

**A wooden mallet to allow for ease of assembly of our interlocking animal models.**

This mallet is made from the same plywood as our models, so won't damage the components during assembly.





## Accessories

### AI Insemination Rod

**SKU**  
**Z1000-000**

This AI training rod is designed for safe and effective practice on Holsim animal models. Made from durable acetyl, it features a 3D printed ergonomic handle for comfort and control, and a soft 3D printed resin tip to prevent puncturing or damaging the silicone during training. Ideal for developing and refining artificial insemination techniques in a realistic, low-risk environment.

## Accessories

### Holsim Universal Stand

**SKU**  
**Z1100-000**

The Holsim Universal Stand is designed to securely hold and display our range of realistic animal models for demonstrations, education, and training. Built with durability and versatility in mind, the stand features dual suction bases with locking levers for a stable, non-slip grip on smooth surfaces. Its modular mounting head allows for easy attachment and positioning of different Holsim animal models, ensuring optimal visibility and access during hands-on use. Made from high-strength materials, this stand is easy to clean, lightweight, and built to withstand frequent transport and setup.

Whether used in veterinary training, classroom demonstrations, or field education, the Holsim Universal Stand provides a reliable and professional platform for interactive learning.



